

ERG 910 PLUS ERG 911 PLUS

Quiet, compact, configurable




SCHILLER
The Art of Diagnostics

ERG 910 PLUS

ERG 911 PLUS

These ergometers perfectly complement SCHILLER exercise ECG and ergo-spiro systems:

- ∴ easy operation via touch display
- ∴ ergonomic features
- ∴ quiet and comfortable
- ∴ compact and stable

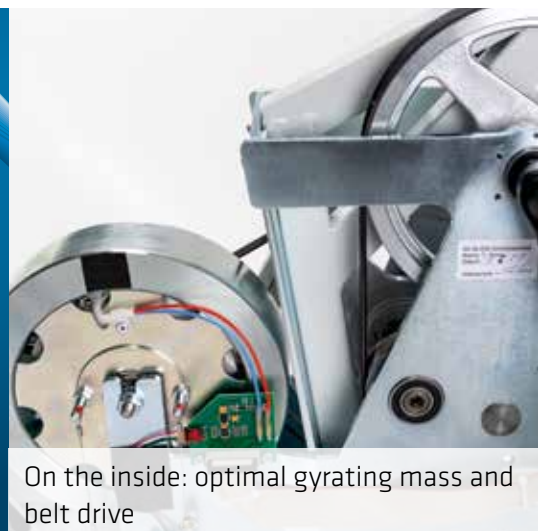

On the inside: optimal gyrating mass and belt drive

The ergometers ERG 910 plus and ERG 911 plus have been developed for use in the cardiology sector. They are very stable thanks to the steel tube construction. Therefore, they are suitable for patients weighing up to 160 kg, with stabiliser (optional) even up to 200 kg. The units are compact and use minimal floor space.

TWO OPERATIONAL MODES

- ∴ Remote operation with a master device (ECG, PC, ...)
- ∴ Standalone operation with either a user-definable step protocol or manual load control

The optimal gyrating mass ensures a comfortable pedalling sensation. Thanks to belt drive, the units are maintenance-free and operate quietly.

ERGONOMIC FEATURES

- ∴ Low frame for easy access
- ∴ Ergonomic sitting position for patients measuring between 1.40 m and 2.10 m, with child saddle and pedals (optional) from 1.20 m
- ∴ Adjustable handle bar (360°)
- ∴ Electrical saddle height adjustment (optional, only for ERG 911 plus)

CONTROL PANEL WITH TOUCHSCREEN

All important exercise data (such as blood pressure, heart rate, pedalling rate, colour-coded information concerning correct pedalling rate) is shown on the display. Blood pressure measurement and electrical saddle height adjustment

are controlled via touchscreen. The control panel is rotatable by 180°.

DIAGNOSTIC AND MONITORING FEATURES

- ∴ SpO₂ measurement with finger or ear clip (optional)
- ∴ Blood pressure measurement with QRS trigger (optional)

CONNECTIVITY

The ergometers are connected to the master device via Bluetooth or RS-232.

ERG 910 PLUS

- ❖ Load range 20-500 watts, optional 20-999 watts

ERG 911 PLUS

- ❖ Load range 20-999 watts
- ❖ Handle bar height adjustable
- ❖ Optional: Saddle height electrically adjustable





Americas
SCHILLER Americas Inc.
 Doral, Florida 33172

North America:
 Phone +1 786 845 06 20
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
 Phone +1 305 591 11 21
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



Asia
SCHILLER Asia-Pacific / Malaysia
 52200 Kuala Lumpur, Malaysia
 Phone +603 6272 3033
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
 A-4040 Linz
 Phone +43 732 70 99 0
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Technology Co., Ltd.
 100102 Beijing, China
 Phone +86 10 84565453
info@schillerchina.com
www.schillermedical.cn



Croatia
SCHILLER d.o.o.
 10000 Zagreb
 Phone +385 1 309 66 59
info@schillerzg.hr
www.schiller.hr



France
SCHILLER Médical
 F-67160 Wissembourg
 Phone +33 3 88 63 36 00
info@schiller.fr
www.schiller.fr



France (distribution France)
SCHILLER France S.A.S.
 F-77608 Bussy St Georges
 Phone +33 1 64 66 50 00
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
 D-85622 Feldkirchen b. München
 Phone +49 89 62 99 81 0
info@schillermed.de
www.schillermed.de



India
SCHILLER Healthcare India Pvt. Ltd.
 Mumbai – 400 059, India
 Phone +91 22 6152 3333 / 2920 9141
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
 PL-02-729 Warszawa
 Phone +48 22 843 20 89 / +48 22 647 35 90
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
 119049 Moscow, Russia
 Phone +7 (495) 970 11 33
mail@schiller.ru
www.schiller.ru



Serbia
SCHILLER d.o.o.
 11010 Beograd
 Phone +381 11 39 79 508
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
 2310 Slovenska Bistrica
 Phone +386 2 843 00 50
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
 ES-28232 Las Rozas/Madrid
 Phone +34 91 713 01 76
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER-Reomed AG
 CH-8912 Obfelden
 Phone +41 44 744 30 00
sales@schiller-reomed.ch
www.schiller-reomed.ch



Turkey
SCHILLER TÜRKİYE
 Okmeydanı-Sisli – Istanbul
 Phone +90 212 210 8681 (pbx)
 Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK Ltd.
 Anstruther, Fife KY10 3H
 Phone +44 16 190 72 565
sales@schilleruk.com
www.schilleruk.com

SCHILLER reserves the right to make changes in design, specifications and features shown, or discontinue the product described at any time without notice or obligation. Images may be representative. Some products may not be regulatory cleared or released for sale in all markets. Contact your SCHILLER representative for the most current information and availability.