

FRED PA-1

Simplicity as a priority



FRED PA-1

FEATURES


- 200J fully automatic defibrillator
- Immediate Power On
- Immediately advises to call emergency number written on the device
- Easy 3 step operation: Start -> Analysis -> Shock
- Highly efficient Lithium batteries deliver 140 shocks at maximum energy
- Automatic energy level selection, required for efficient shock delivery
- SD Card facility to store data during intervention
- Capacity of storing 500 events and 2hrs of ECG recording storage in the SD card
- Complete user guidance by voice prompts & pictograms for easy operation
- Metronome facility for accurate CPR
- LED indications for Maintenance, Replace Electrode and Replace Battery

EASY ACCESS

Electrode holding slot

Call

Emergency number

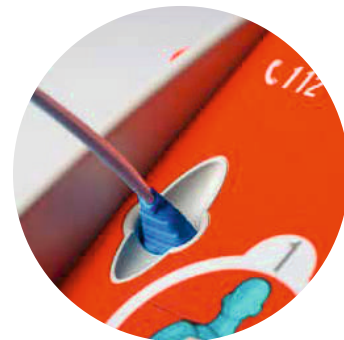


No time delay

Pre-connected electrode allows faster analysis

Safe & Ready to use

Self-test indicator with ready-to-use indication



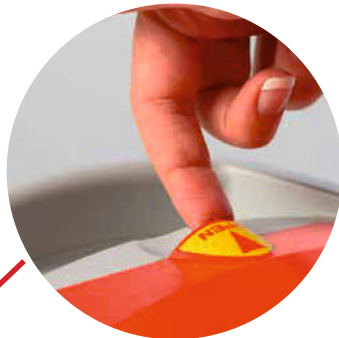

Saves time

Auto-switch ON once the cover is lifted.
Time efficient device fit to save a life



Repeat

Repeat last instruction given to user
in case an instruction is missed out

Warnings

LEDs indicating status of the device



User friendly

Easy user interface with Pictograms
give clear 1-2-3 step guide



SCHILLER, founded in 1974 by Alfred E. Schiller is a successful group with 30 subsidiaries and a global sales network. Today, SCHILLER is a world-leading manufacturer and supplier of devices for cardiopulmonary diagnostics, defibrillation and patient monitoring as well as software solutions for the medical industry.

Over the last 27 years, SCHILLER has established itself as a leader in the medical technology sector in India. With 400+ employees, two production centres and a R&D facility SCHILLER makes state-of-the-art healthcare equipment accessible through a network of 100 sales & service personnel spread across more than 45 locations. Our product range covers Critical Care, Anaesthesia, Emergency Care, Cardiology, Pulmonology, Radiology and Robotics.

Scan here
to know more



Part No. B4.100007 9/2024 Ver. 1

All registered trademarks acknowledged.

SCHILLER
The Art of Diagnostics

Swiss H.Q.: SCHILLER AG, Altgasse 68, P. O. Box 1052, CH -6341 Baar, Switzerland

Indian H.O.: Schiller Healthcare India Pvt. Ltd., Advance House, Makwana Road, Andheri (East), Mumbai-59. Tel.: + 91 - 22 - 61523333

Factory: No. 15/5 & 15/6, Vazhuthavur Road, Kurumbapet, Puducherry - 605 009, India

Regional Offices: Mumbai: 022-61523333 Email: sales.west@schillerindia.com Delhi: 011-41062067 / 09312432205 Email: sales.north@schillerindia.com Kolkata: 033-23593102 / 033-23593103 Email: sales.east@schillerindia.com Chennai: 044-28232648 / 044-28311021 Email: sales.south@schillerindia.com Ahmedabad: 079- 35337949 / 079- 35337950 Email: sales.west@schillerindia.com Puducherry: 0413- 220 2680 / 0413- 229 2940 Email: sales.south@schillerindia.com

CIN : U33110MH1997PTC111307

sales@schillerindia.com, help@schillerindia.com www.schillerindia.com 1800 209 8998

Follow us on : /SchillerHealthcareIndia /schiller_healthcare_india /schillerhealthcareindia /schiller_healthcare_india @schillerindia1074